

⚠ Caution

- Do not use water to wash the unit.
This may cause error



- This unit is only for measure weight.
Do not use for any other purpose



- If user act inconvenience when use this unit,
please ask other help or use armrest.




- The unit can't measure correctly on carpet.
Make sure you put on flat floor to use.




Indication

- Use the unit in assign battery type.
Do not mix old and new one together.

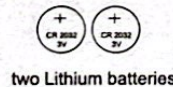

- When install battery, do not place wrong
the +(positive and - (negative) direction

- Remove the batteries if you do not intend to
use the unit for a period of three months or
longer before placing the device in storage

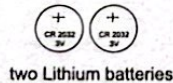

- Dispose of used batteries according to applicable
local regulations.

chorder HEBEI

Product Feature

- Highly precision: minimum display unit is 0.1kg. The slight weight variation can be measure.
- Tap On SW function: only step on softly the unit left below or light below can turn on.
- Calculate and display user BMI.
- Storage four user's height setting and also storage last time measure result.
- Large LCD display.
- Automatically turn off: approximately six seconds automatically turn off after measure.

Know how to use correctly

■ Install battery

Install two 3V lithium batteries before use the unit
This unit have attached trial battery,
if the battery is low battery, you need to change

- Open the battery cover
- Install four 1.5V dry batteries into battery box (differentiate between + (positive and - (negative)), replace the battery cover

■ How to take a measurement

- Step on softly the unit left below or light below. After turn on, your foot have to leave the unit (if your foot keep step on the unit, it will not be measure correctly)
- Display " _ _ _ _ " Wait for the 0.0 symbol to appear before step on the unit for measure.(if the 0.0 symbol to appear before you step on the unit, it will not measure correctly)
- Start proceeding, after display " _ _ _ _ " shows measure result. The monitor automatically turn off after six second. if after " _ _ _ _ " display, the measuring is not stable, it will display " _ _ _ _ " until the measure complete, shows measure result. If circulate display over 10 second, it will not shows measure result. It display "Err" and then turn off after three seconds.

■ Setting Personal Data

1	Press SET button. The user code U-1 blinks on the display
2	Select the UP button or Down to adjust personal profile number (1~4), or Guest mode (user9). The Guest mode don't have memory function, only for measure.
3	Press the SET button to advance height setting. The default value "170" blinks on the display. Select the UP button or Down to change the height. Press and hold the button to advance at a faster speed
4	Press the SET button to advance measure mode. (refer to measurement 2,3 step). In the Guest mode, after display shows ZERO, you can Select the UP button or Down to check last time data
5	The display automatically repeatedly three times (1.5 per second) to shows Weight and BMI ,and then automatically turn off after nine second
6	When display automatically repeatedly Weight and BMI, you can Select the UP button or Down to check the last time Weight and BMI

■ Deleting Personal Profile

Turn the monitor off. Press and hold the UP button and Down and SET button, the monitor automatically turn on and display CLR symbol. This means the USER1~USER4 personal profile is deleted from the memory

ERROR INDICATORS

ERROR DISPLAY	CAUSE	CORRECTION
E	weight exceeds 150 kg	
Err	you stepped onto the unit all along after turn on	do not stepped onto the unit all along after turn on

Battery Information

- (1) When the LO Battery Indicator appears on the display screen, replace with new battery.
- (2) Use two 3V lithium Battery.
- (3) Old batteries may not even display the LO battery indicator.
- (4) Battery life approximately six months (when used three times a day).
- (5) Remove the battery to storage if you do not intend to use the unit for a period . Besides install the battery not use the unit for a long period, the battery may leak liquid and cause trouble.

Care and Maintenance

Avoid using any gas oil, lacquer thinner, hot water to clean.
Clean the unit dirt with a soft cloth soak in a neutral clean liquid,
and then use other soft cloth to clean up

If you put article on the unit for a long time, the unit will be power on condition.
This cause battery consumption faster result.

Troubleshooting Tips

Problem	Cause	Solution
The unit can not turn on	the unit is place the 凹凸 floor to measure	make sure measure on the flat floor
	the battery is consumption	change to new battery
	the battery +(positive and - (negative) is not in correct direction	place right direction
he weight value is displayed abnormally low or high	the unit is place mat or carpet to measure	make sure measure on the flat floor
	the unit is place the 凹凸 floor to measure	do not stepped onto the unit all along after turn on
Other abnormally display	SW make wrong operation	remove the batteries and reinstall

Specifications

Measure range	5~150 kg	Battery life	approximately six months (when used three times a day)
Minimum display unit	0.1 Kg	External dimension	approximately 316 X approximately wide271 X approximately 41mm
Power	2X3V Lithium Button Battery(2 PCS of Trial Battery attached)	Weight	approximately 1.2kg (not included batteries)
Attachment	Instruction manual, two trial battery.		

What is the BMI?

BMI is "Body Mass Index.(weight in kilogram) / (height in metre) / (height in metre)